HYDRA - HIGH PURITY

HYDRA SERIES

HIGH PURITY
GENERAL DEIONIZATION
INDUSTRIAL & LABORATORY

The Hydra High Purity DI System is a pre-assembled cartridge system that provides deionized water using a staged filtration process. Water purification is provided using a 3-step process. Stage (1) removes particles greater than 5 micron, chlorine, and organics. Stages (2 & 3) are deionizers designed to remove dissolved minerals. The Hydra High Purity DI series includes all cartridges and resistivity light for a total turnkey system. The Hydra high purity DI configuration uses low oder, Type II anion mixed bed cartridges.



HIGHLIGHTS

- Wall Mounted and Low Profile
- Water Quality Light Indicating Cartridge Replacement
- Economically Produces DI Water
- Turnkey System For On Demand Service
- Designed For City Water & Tap Water Supply
- Outlet Isolation Valve
- Spanner Wrench & Associated Tubing Provided

SYSTEM SPECIFICATIONS

Pressure: 20 - 90 PSIG (changed values)

Temperature: 100° F Max

Voltage: 120 V Connections:

> Inlet – 3/8" OD Tubing Outlet – ½" Hose Barb

Materials:

Housing – Polypropylene Bracket – Painted Steel

APPICATIONS

- General deionization
- Battery water filling
- Humidification
- Hydrogen generator
- Glassware rinse
- Glassware washer
- Sterilizers





HYDRA DI CAPACITY

| FEEDWATER TOTAL DISSOLVED SOLIDS (TDS AS CaCO₃) | HYS-122-S-HP (Gallons) | HYS-122-B-HP (Gallons) |
|---|---------------------------|---------------------------|
| 10 ppm* | 2,650 | 6,800 |
| 300 ppm | 90 | 225 |


^{*} Typical water quality with reverse osmosis pretreatment.

ORDERING GUIDE

| BOWL DIAMETER | SHIPPING WEIGHT | PART NUMBER (Legacy Part Number) | REPLACEMENT KIT PART NUMBER (Legacy Part Number) | REPLACEMENT RESISTIVITY LIGHT 1 MΩ-cm |
|------------------|--------------------|-------------------------------------|--|---|
| 2.5" | 24 lbs. | HYS-122-S-HP (HY-122-DI) | HYK-122-S-HP (HYK-001) | CL-200K-50 |
| 4.5" | 55 lbs. | HYS-122-B-HP (HY-122B-DI) | HYK-122-B-HP (HYK-002) | CL-200K-50 |

SPECIFICATIONS

| PART NUMBER (as Listed on Website) | DIMENSIONS (H X W X D) | FLOW RATE | TOTAL CAPACITY GRAINS AS CaCO₃ | TYPICAL EFFLUENT |
|------------------------------------|---------------------------|-----------|-----------------------------------|---------------------|
| HYS-122-S-HP | 25"x 20"x 6" | 0.5 gpm | 1,800 | 10 MΩ-cm |
| HYS-122-B-HP | 26"x 29"x 9" | 1.25 gpm | 4,760 | 10 MΩ-cm |

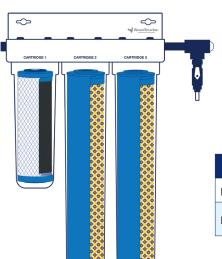

| | CARTRIDGE 1 | CARTRIDGE 2 | CARTRIDGE 3 |
|-------------|--|--|--|
| DESCRIPTION | Carbon Pretreatment | Deionization | Deionization |
| DIMENSIONS | 2.5" X 10" (Slim) 4.5" X 10" (Big Blue) | 2.5" X 20" (Slim) 4.5" X 20" (Big Blue) | 2.5" X 20" (Slim) 4.5" X 20" (Big Blue) |

Revision 1.0 © 2020 ResinTech, Inc.

